

# JAY<sup>®</sup>

# ION CUSHION

## JAY<sup>®</sup> Ion<sup>™</sup> Cushion

### IMPORTANT CONSUMER INFORMATION

**NOTICE:** This manual contains important instructions that must be passed on to the user of this product. Please do not remove this manual before delivery to the end user.

**SUPPLIER:** This manual must be given to the user of this product.

**USER:** Before using this product, read this entire manual and save for future reference.

Owner's Manual

## Cojín JAY<sup>®</sup> Ion<sup>™</sup>

### INFORMACIÓN IMPORTANTE PARA EL CONSUMIDOR

**AVISO:** Este manual contiene instrucciones importantes que deben ser entregadas al usuario de este producto. Por favor no retire este manual antes de la entrega al usuario.

**DISTRIBUIDOR:** Este manual debe ser entregado al usuario de este producto.

**USUARIO:** Antes de usar este producto, lea este manual en su totalidad y guárdelo para futura referencia.

Manual de instrucciones

## Coussin JAY<sup>®</sup> Ion<sup>™</sup>

### INFORMATIONS IMPORTANTES DESTINÉES AU CONSOMMATEUR

**AVIS :** Ce manuel contient des instructions importantes qui doivent être communiquées à l'utilisateur de ce produit. Prière de ne pas le retirer avant livraison.

**FOURNISSEUR :** Ce manuel doit être remis à l'utilisateur de ce produit.

**UTILISATEUR :** Avant d'utiliser ce produit, lisez entièrement ce manuel et conservez-le pour référence ultérieure.

Manuel d'instructions



MACK'NME



**NOTE:** Check all parts for shipping damage. In case of damage DO NOT use. Contact Carrier/Sunrise for further instructions.

Sunrise Medical recommends that a clinician such as a doctor or therapist experienced in seating and positioning be consulted to determine if a JAY Ion Cushion is appropriate. Cushions should only be installed by an authorized Sunrise Medical supplier.

## GENERAL WARNINGS

### WARNING

**DO NOT** install this equipment without first reading and understanding this manual. If you are unable to understand the Warnings, Cautions and Instructions, contact a qualified clinician or supplier - otherwise injury or damage may occur.

### WARNING

Hook and Loop velcro should not come into prolonged contact with moisture as this may degrade the adhesive and lead to a failure of the velcro.

### WARNING

Installing a cushion on a wheelchair may affect the center of gravity of the wheelchair and may cause the wheelchair to tip backwards, potentially resulting in injury. Always assess for the potential need for antitippers or amputee axle adapter brackets to be added to the wheelchair to help increase stability.

### WARNING

Prior to prolonged sitting, any cushion should be tried for a few hours at a time while a clinician inspects your skin to ensure that red pressure spots are not developing. You should regularly check for skin redness. The clinical indicator for tissue breakdown is skin redness. If your skin develops redness, discontinue the use of the cushion immediately and see your doctor or therapist.

The JAY Ion cushion is designed to help reduce pressure. However, no cushion can completely eliminate sitting pressure or prevent pressure sores. The JAY Ion cushion is not a substitute for good skin care including: proper diet, cleanliness, and regular pressure relief.

**CAUTION-** AVOID SHARP OBJECTS OR EXPOSURE TO EXCESSIVE HEAT OR OPEN FLAME.

### WARNING

**OBSTRUCTIONS: DO NOT** place any obstructions between the user and the cushion because it will reduce product effectiveness.

### WARNING

**COVER ORIENTATION:** Product must be used with the cushion cover facing up. If the cover is not used correctly it may reduce or eliminate the cushion's benefits and could increase risk to the skin and soft tissue.

## MAXIMUM USER WEIGHT:

14" to 21" width = 300 lbs.

22" to 24" width = 500 lbs.


## JAY ION CUSHION

The JAY Ion is designed with the benefits of skin protection and positioning in a lightweight, comfortable and low-maintenance design. The JAY Ion incorporates an antimicrobial outer cover, a moisture resistant inner cover and a dual layer foam base.

The JAY Ion helps to provide pressure and shear reduction which may help reduce the incidence of decubitus ulcers. The cushion is intended to be comfortable and does not restrict movement, allowing the user to change positions on the cushion.

The maximum user weight for this product is 300 lbs. for 14" to 21" widths and 500 lbs. for 22" to 24" widths.

## JAY ION COMPONENTS

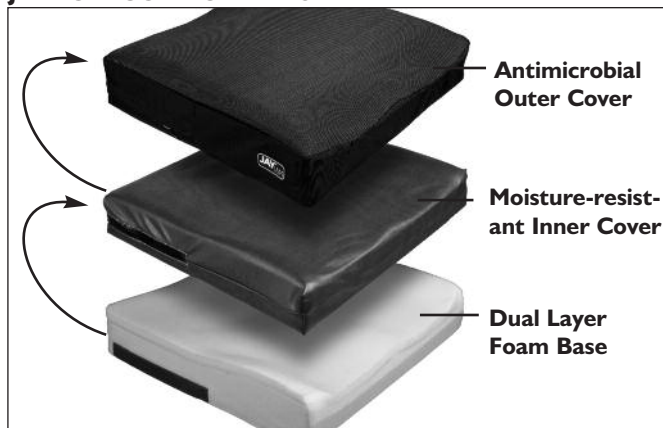

## FITTING A PRESSURE DISTRIBUTING CUSHION

Select a cushion to match your hip width when in the seated position. The cushion length should extend to within 1" - 2" (2.5cm - 5cm) from the back of the knees. Place the cushion to the back of the wheelchair with the pelvic loading area toward the rear. When fit correctly, your hips should be to the back of the chair. Your ischials (seat bones) should be centered in the rear of the cushion. Proper footrest adjustment may enhance sitting comfort and help lower peak sitting pressures. Adjust the footrests so your legs rest firmly but comfortably on the cushion.

Avoid using the cushion on overly stretched seat upholstery. Stretched upholstery does not provide adequate support for the cushion and improper positioning may result. Replace stretched upholstery or purchase the JAY Adjustable Solid Seat.

Avoid using overly stretched backrest upholstery. Stretched back upholstery does not provide adequate support for the body and improper positioning may result. Replace stretched backrest upholstery or purchase a back support product from Sunrise that is appropriate for you. Contact an authorized supplier for a free demonstration of these products.

## EASY MAINTENANCE & CLEANING

Regular cleaning and maintenance may help extend the life of your cushion. During cleaning, component inspection is recommended.

Check the cover for tears and excessive wear and or any other abnormalities. Check the foam base to ensure foam consistency.

### To clean cover

1. Remove the outer cover from the foam base and zip closed.
2. Machine wash in warm water and drip or tumble dry on low heat.

No ironing required.

**NOTE:** Do not dry clean covers or use industrial washers and dryers to clean cover. Do not steam autoclave. Do not bleach.

To refit cover after washing, fit cover back over foam base. Ensure the back of the cover matches up to the back of the foam base.

**NOTE:** If utilized in an institutional setting, write the resident's name with a permanent marker on the content label of the covers. This will help to ensure the return of the covers after washing.

### To clean foam base

Remove the cover from the foam base. Wipe lightly with a damp cloth. Do not use soap. Do not submerge in water. Wipe off with a clean cloth. Let it completely air-dry before reassembly.

**NOTE:** Avoid submersion of the foam base in water.

## REASSEMBLING THE JAY ION CUSHION AFTER CLEANING

Follow cleaning and maintenance instructions as described in "Easy maintenance and Cleaning" then reassemble as follows.

Insert the foam base into the inner cover with white label down and zipper to rear of cushion. Then insert base plus inner cover into the outer cover. Make sure white inner cover label is down and zipper is to the rear. Ensure that both zippers are at the rear of the foam base (the part number tag should be to the rear of the foam base) where the pelvis is typically located.

## JAY ION CUSHION WARRANTY

Every JAY Ion cushion base is guaranteed to be free from defects in materials and workmanship for a period of 24 months from the date of purchase, provided normal use. Should a defect in materials or workmanship occur within 24 months from the original date of purchase, Sunrise Medical will, at its option, repair or replace it without charge. This warranty does not apply to punctures, tears or burns, nor to the cushion's removable cover.

The removable outer and inner covers are guaranteed to be free from defects in materials and workmanship for a period of six (6) months from the date of purchase, provided normal use. Should a defect in materials or workmanship occur within six (6) months from the original date of purchase, Sunrise Medical will, at its option, repair or replace it without charge.

Claims and repairs should be processed through the nearest authorized supplier. Except for express warranties made herein, all other warranties, including implied warranties of merchantability and warranties of fitness for a particular purpose are excluded.

There are not warranties which extend beyond the description on the face hereof. Remedies for breach of express warranties herein are limited to repair or replacement of the goods. In no event shall damages for breach of any warranty include any consequential damages or exceed the cost of non-conforming goods sold.



# MACK'NME

## Mobility Solutions

144 Lake Entrance Rd  
Oak Flats, NSW, 2529  
P:02 4257 1976  
M:0417 106 370  
W:macknme.com.au  
E:admin@macknme.com.au

### Sunrise Medical (US) LLC

2842 Business Park Ave.  
Fresno, CA 93727  
USA  
(800) 333-4000  
In Canada (800) 263-3390  
www.sunrisemedical.com

### Sunrise Medical, Inc.

237 Romina Drive, Unit 3  
Concorde, Ontario L4K 4V3  
CANADA  
(800) 263-3390  
www.sunrisemedical.ca

### Sunrise Medical Limited

High Street Wollaston  
West Midlands DY8 4PS  
ENGLAND  
Phone: +44 (0) 13 84 44 66 88  
www.sunrisemedical.com

### Sunrise Medical GmbH+Co.KG

69254 Malsch/Heidelberg  
DEUTSCHLAND  
Tel.: +49 (0) 72 53/980-0  
Fax: +49 (0) 72 53/980-111  
www.sunrisemedical.com

### Sunrise Medical S.L.

Polígono Bakiola, 41  
48498 Arrankudiaga – Vizcaya  
ESPAÑA  
Tel.: +34 (0) 9 0214 24 34  
Fax: +34 (0) 9 46 481575  
www.sunrisemedical.es

### Sunrise Medical S.A.

13, Rue de la Painguetterie  
37390 Chanceaux sur Choisille  
FRANCE  
Tel.: +33 (0) 247 55 44 00  
Fax: +33 (0) 247 55 44 03  
www.sunrisemedical.com

### Sunrise Medical Srl

Via Riva 20, Montale  
29100 Piacenza (PC)  
ITALY  
Tel.: +39 (0) 05 23-57 3111  
Fax: +39 (0) 05 23-570 60  
www.sunrisemedical.it

### Sunrise Medical A.G.

Lückhalde 14  
3074 Muri bei Bern  
SCHWEIZ  
Tel.: +41 (0) 31-958-38 38  
Fax: +41 (0) 31-958-38 48  
www.sunrisemedical.ch

### Sunrise Medical AS

Rehabsenteret  
1450 Nesoddtangen  
NORWAY  
Tel.: +47 66 96 38 00  
Faks: +47 66 96 38 80  
www.sunrisemedical.no

### Sunrise Medical AB

Britta Sahlgrens gata 8A  
421 31 Västra Frölunda  
SWEDEN  
Tel: +46 (0) 31-748 37 00  
Fax: +46 (0) 31-748 37 37  
www.sunrisemedical.se

### Sunrise Medical B.V.

Groningenhaven 18-20  
3433 PE Nieuwegein  
The NETHERLANDS  
Tel.: +31 (0)30 – 60 82 100  
Fax: +31 (0) 30 60-5 58 80  
www.sunrisemedical.com



This product is manufactured to comply with the 'Medical Device Directive' 93/42/EEC.

Questo prodotto è stato fabbricato in conformità alla direttiva 93/42/CEE relativa ai dispositivi medici.

Denna produkt har tillverkats i enlighet med EU:s direktiv om medicintekniska produkter nummer 93/42/EEC.

Dieses Produkt ist in Konformität mit der 'Richtlinie für Medizinprodukte' 93/42/EEC hergestellt.